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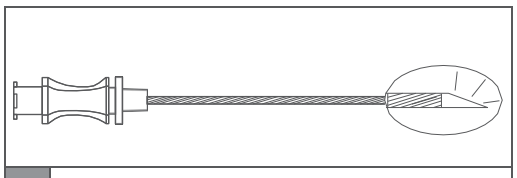
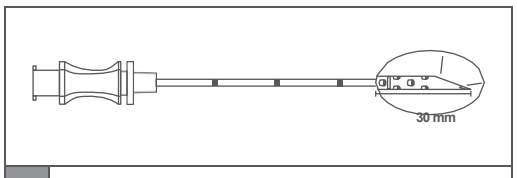
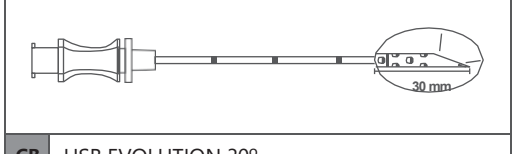
GB

Instructions for use

Only trained and experienced health care professionals should use this equipment. Before using any system component or any component compatible with this system, read and understand the instructions. Pay special attention to WARNING information. Failuretocomplymayresultininjuriesofpatientand/or operating room staff. Upon initial receipt and before each use, examine system components and inspect device for damage. Do not use if evident damage.

INDICATIONS

Device is designed for peripheral nerve location for regional anaesthesia with the help of an ultrasound medical device only.

CONTRAINDICATIONS

- Not intended for neurolytic ablation.
- Not intended for epidural use.
- Do not use if any component is damaged or missing or pac-kaging damaged.

WARNINGS

- In the case of bone contact, the needle must be withdra replaced by anyonewn and examined, if it has been damaged it must be replaced by a new one.
- Needle is for use with a low frequency ultrasound devices only.
- Sterile only if package is unopened and undamaged.
- For single use only.
- DO NOT reuse or resterilize.

The health care provider performing any procedure is responsi-ble for determining the appropriateness of this instrument and the specific technique for each patient. Temena, as manufacturer, DOES NOT recommend surgical procedure.

COMPLICATIONS

Though rare, nerve blocks have been known to cause some side effects and can sometimes result in serious health complications. Possible side effects include:

- Paralysis
- Paresthesias
- Nerve damages
- Cardiac arrest (very rare)
- Neurologic complications
- Pain around the injection site
- Inadvertent puncture of the arteria

SAFETY MEASURES

- Before insertion, correct the patient's hemodynamic deficiency.
- Use aseptic techniques
- Open well-dimensioned venous accesses (16 G -18 G) in order to compensate hemodynamic changes.
- Have vasoconstrictors prepared and available (e.g. Ephedrine)
- Have all material for intubation and artificial resuscitation as well as for cardiopulmonary resuscitation available.

INSTRUCTIONS FOR USE

1. Use aseptic techniques.
2. Prep and cover the patient with a sterile drape in accordance with standard institutional procedures for the block techni-que.
3. Cover the ultrasound probe within a sterile bag, put a sterile gel on the skin according to the instruction for use.of the ultrasound device's
4. To obtain the best visibility of the needle the ultrasound pro-be should be oriented in a transversal manner relative to the needle. This means the ultra sound beam should hit the target in an angle of approximately 90° or the targeted nerve.
5. Insert the needle.
6. Check if the needle's tip and the echogenic area of the needle is visible on the display.of the ultrasound device's.
7. Advance the needle until it can be seen close to the target nerve (s). Note: You can inject between 0.5 and 2 ml of liquid by bolus or by continuous infusion, this way creating a small hypoechogenic area that may help to line up the needle and the probe.
8. Connect the syringe to the needle and after an aspiration test inject cold physiological solution generating paresthesias
9. Inject one to three millilitres of local anaesthetic solution by intermittent injection. Reposition the needle at least once under the ultrasound observation. The elicited motor respon-se should extinguish within 10 seconds. (Continued motor response may indicate intravascular or intramuscular needle placement.)
10. Inject progressively the remaining local anaesthetic solution, aspirating frequently.
11. After completion of this procedure, dispose all components in the appropriate container.

LIMITATION

In some cases:
Nerves can move to the side when approached by the needle. Fascial membranes could drain out the local anaesthetic from the target nerve. To elicit muscle contraction, it is necessary to make a minor adjustment by moving the needle longitudinally along the course of the nerve.

NOTE

TEMENA NEEDLES HAVE BEEN DESIGNED TO PROVIDE THE VERY BEST RESULTS. USE CAREFULLY AND ACCORDING TO THE INSTRUCTIONS ABOVE. ALL US NEEDLES HAVE BEEN DESIGNEDFORSINGLEUSEANDHAVEBEENSTERILIZEDBYETHYLENE OXYDE. IT IS NOT ALLOWED TO RESTERILIZE THE NEEDLE, BECAUSE THE PRODUCT CAN BE DESTROYED. THIS KIT IS PHTHALATES- FREE. DO NOT USE IF THE PACKAGING HAS BEEN DAMAGED.

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